


Dräger DrugCheck® 3000 Drug Testing Device

Use the Dräger DrugCheck® 3000 to find out within minutes if a person recently consumed certain drugs. The compact and quick saliva (oral fluid) based drug test yields reliable results affordably and easily. The non-invasive drug test kit does not require electricity and can therefore be used anywhere.

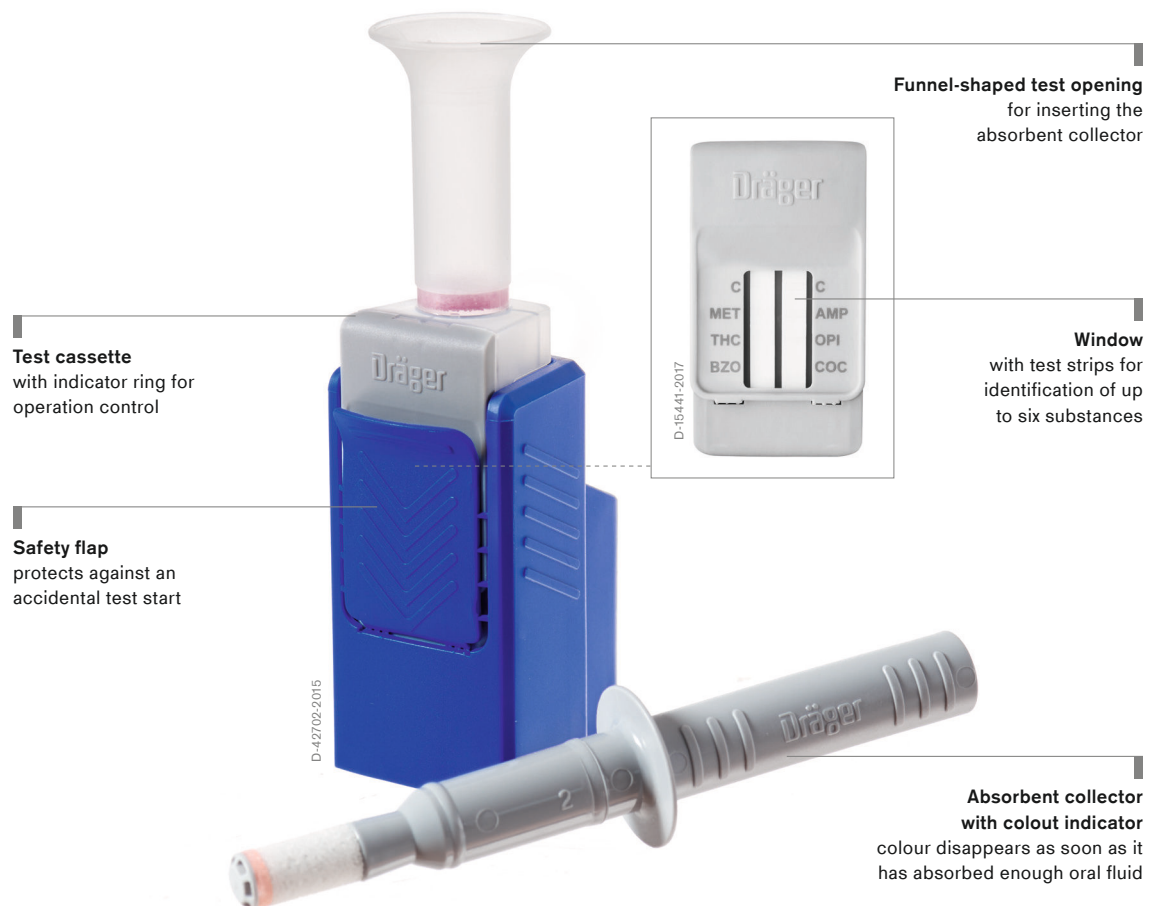

Benefits

Drug test – obtain samples easily and safely

The Dräger DrugCheck 3000 consists of two components: an absorbent collector for obtaining a saliva sample, and a test cassette for the analysis. The test cassette contains the buffer liquid and a window with two test strips, which displays the control and test lines. The drug test itself is performed in three easy steps: swab the saliva sample, shake the test kit, wait for a brief incubation period, and then start the test. A colour indicator on the swab disappears as soon as it has absorbed enough oral fluid for a test.

As soon as the control lines appear in the window, you can read the results. If the test result is negative, a line will appear alongside the respective substance class (drug). This means that none of the target substance was detected in the sample. If a line fails to appear next to one of the substance classes, then the result for this substance is positive.

Fast and sensitive drug analysis

Check individuals for up to six substance classes simultaneously with Dräger DrugCheck 3000: cocaine, opiates, amphetamine, methamphetamines/designer drugs (e.g. Ecstasy, MDMA), benzodiazepines (e.g. in prescription medicines) and cannabis (THC). Of all the substance classes listed, cannabis is the drug consumed most frequently and also the most difficult of all compounds to identify (THC = delta-9-tetrahydrocannabinol). This is why the Dräger DrugCheck 3000 was optimised to detect THC, and now offers two measurement options: depending on the chosen waiting time (pre-incubation) of the test there is a fast option with a higher or sensitive option with a low THC cut-off. The entire drug test only takes 3 – 5 minutes.

Unambiguous on-the-spot drug screening

Its compact pocket sized design makes the DrugCheck 3000 easy to transport. It can be made ready for use quickly and easily, and is suitable for flexible use applications such as roadside checks and testing in safety-relevant workplaces*. The test kit has no electrical parts, which makes the drug tester safe to use in hazardous areas and even as point-of-care-testing (POCT).

Additional benefits

Dräger has decades of experience in alcohol measuring equipment and drug detection methods. For the DrugCheck 3000, Dräger employed the testing principle of the proven Dräger DrugTest® 5000 system, which reliably detects even minute traces of THC. The disposable DrugCheck 3000 test kit cannot be manipulated and is hygienic to use.

*The DrugCheck 3000 is sold only for law enforcement purposes in the USA.

Related Products



D-54720-2012

Dräger DrugTest® 5000

No pipetting, no drips, no timing: With the Dräger DrugTest® 5000, a drug test is carried out simply and quickly. The collected saliva sample can be analysed immediately - for accurate results on the spot.



D-6296-2009

Dräger SSK 5000

Ideal for surface drug testing applications, the Dräger SSK 5000 sampling system is quick and easy to use. When used with the Dräger DrugTest® 5000, the sample can be transferred to a test cassette and screened for possible drug contamination. The sample can also be stored and transported for further testing.



ST-13226-2007

Dräger DCD 5000

Simple, fast and secure detection of intoxicants and drugs: The Dräger DCD 5000 does not only enable an easy collection of saliva, you can also store the samples and send them to the laboratory for analysis – without the sample being falsified along the way.



D-4347-2014

Dräger SCK

To follow the chain of custody please use the confirmation kit with oral fluid collector for subsequent sampling for laboratory analysis of one sample including shipping material.

Technical Data

Test for detecting drugs in oral fluid

Dimensions (W x H x D)	32 x 111 x 57 mm 1.26 x 4.37 x 2.24 inches
Weight	< 30 g (0,07 lb)
Operating range	Operation: 5 °C to 40 °C (41 °F to 104 °F) at 5 % to 95 % relative humidity
Storage / transportation	4 °C to 25 °C (39 °F to 77 °F)
Time of measurement	Fast measurement: analysis < 3 minutes Sensitive measurement: analysis < 5 minutes
Selecting the measurement mode	Depending on the selected pre-incubation time and the desired THC cut-off fast or more sensitive reading.
Licenses / standards	Licensed as a medical product within the EU in accordance with Directive 98/79/EC on in vitro diagnostic medical devices (IVD). Outside of the European Union for non-medical use (Non-IVD) or for forensic use (only USA).

Ordering Information

Dräger DrugCheck® 3000	Description	Ordering no.
Test for detecting drugs in oral fluid		
Dräger DrugCheck® 3000 STK 5 IVD	Packaging unit 20 kits	83 25 500
Dräger DrugCheck® 3000 STK 5 Non-IVD	Packaging unit 20 kits	83 25 580
The Dräger DrugCheck® 3000 STK 5 detects the following substances:		
<ul style="list-style-type: none"> - Cocaine (COC) - Opiates (OPI) - Cannabis (THC) - Amphetamines (AMP) - Methamphetamines (MET) 		
Dräger DrugCheck® 3000 STK 6 IVD	Packaging unit 20 kits	83 27 950
Dräger DrugCheck® 3000 STK 6 Non-IVD	Packaging unit 20 kits	83 27 960
The Dräger DrugCheck® 3000 STK 6 detects the following substances:		
<ul style="list-style-type: none"> - Cocaine (COC) - Opiates (OPI) - Cannabis (THC) - Amphetamines (AMP) - Methamphetamines (MET) - Benzodiazepines (BZO) 		

Notes

Notes

Not all products, features, or services are for sale in all countries.
Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
www.draeger.com

REGION DACH

Dräger Safety AG & Co. KGaA
Revalstraße 1
23560 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION EUROPE

Dräger Safety AG & Co. KGaA
Revalstraße 1
23560 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION MIDDLE EAST, AFRICA

Dräger Safety AG & Co. KGaA
Branch Office
P.O. Box 505108
Dubai, United Arab Emirates
Tel +971 4 4294 600
Fax +971 4 4294 699
contactuae@draeger.com

REGION ASIA PACIFIC

Draeger Singapore Pte. Ltd.
25 International Business Park
#04-20/21 German Centre
Singapore 609916
Tel +65 6308 9400
Fax +65 6308 9401
asia.pacific@draeger.com

REGION CENTRAL AND SOUTH AMERICA

Dräger Panama S. de R.L.
Complejo Business Park,
V tower, 10th floor
Panama City
Tel +507 377-9100
Fax +507 377-9130
contactcsa@draeger.com

Locate your Regional Sales
Representative at:
www.draeger.com/contact

